

BiFlex

URETERAL ACCESS SHEATH

More Than An Access Sheath

Designed to ease navigation through a tortuous urinary tract, the BiFlex Ureteral Access Sheath has a dilating catheter with two working channels. At 3 Fr each, the channels include luer lock connections, allowing for a safety wire and the injection of a contrast medium.



Cogentix
Medical

BiFlex

URETERAL ACCESS SHEATH

» 2 Working Channels (3 Fr)

For the use of a safety wire and injection of contrast medium

» Hydromer® Hydrophilic Coating

Facilitate the introduction and the progression through a tortuous urinary tract

» Kink-Resistant Coil

Allows easy entry and withdrawal of the endoscope during procedure

» Distal Platinum Ring

Confirms location of distal tip under fluoroscopy



» Tapered Introducer / Dilator

Tapers to 4.8 Fr for an atraumatic insertion/dilation

» Large Inner Diameter

“... The inner diameter of the 10/12 Bi-Flex is actually 11 F, so intrarenal pressures were lower and irrigation flow was higher ... The irrigation flow measurements were in accordance with pressure measurements.”¹

1. T. Sener, J. Cloutier, L. Villa, F. Marson, Buttice, S. Doizi and O. Traxer
France Journal of Endourology, September 18, 2015

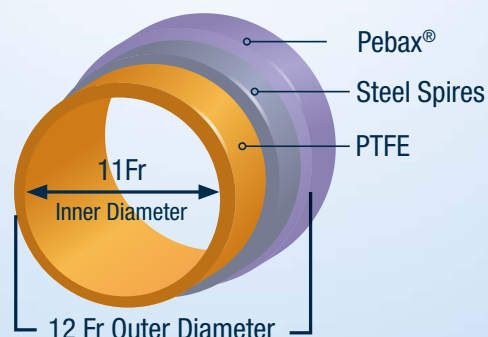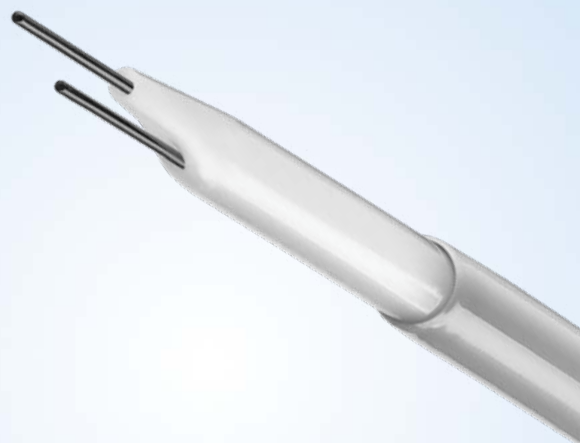

Product Information

CATALOG NUMBER	DESCRIPTION
09-0015	Ureteral access sheath with dual lumen dilator (3 Fr working channels), 10/12 Fr x 35cm
09-0016	Ureteral access sheath with dual lumen dilator (3 Fr working channels), 12/14 Fr x 35cm
09-0017	Ureteral access sheath with dual lumen dilator (3 Fr working channels), 10/12 Fr x 45cm
09-0018	Ureteral access sheath with dual lumen dilator (3 Fr working channels), 12/14 Fr x 45cm

Cogentix
Medical

Cogentix Medical
5420 Feltl Road
Minnetonka, MN 55343

cogentixmedical.com

TEL 866.258.2182
FAX 866.255.4522